



# biomarkers in the **management** of **CARDIAC EMERGENCIES**

- Acute Coronary Syndromes
- Acute Heart Failure





## CARDIOVASCULAR DISEASE, a major public health issue

### ● ACUTE CORONARY SYNDROMES (ACS):

- First cause of death worldwide, <sup>(1)</sup>
- Leading cause of disease burden in high-income countries, <sup>(1)</sup>
- Over 5 million annual hospital discharges in Europe and the USA, <sup>(2, 3)</sup>
- Total annual cost of USD 215 billion in Europe and the USA. <sup>(2, 4)</sup>

### ● HEART FAILURE (HF):

- High proportion (1 in 3) of individuals aged 55 will develop HF during their remaining lifespan, <sup>(5)</sup>
- Total annual cost of USD 100 billion in Europe and the USA, 70% of which is due to hospitalization. <sup>(6)</sup>

**Rapid exclusion** is particularly useful in ED patients with chest pain and/or acute dyspnea, because the majority will not have ACS or acute HF:

- Ischemic cardiac disease is present in less than 25% of patients with suspected ACS, <sup>(7)</sup>
- Acute HF is present in 35% of patients with acute dyspnea. <sup>(8)</sup>

*For list of abbreviations, refer to page 31.*





# Introduction

**Chest pain** and **shortness of breath (SOB)/dyspnea** often occur together and are among the most frequent complaints in patients visiting an emergency department (ED)<sup>(8,9)</sup>.

The evaluation of patients with these symptoms is a **challenge** for the ED physician because of the variety of potential causes and comorbid medical conditions.

It is particularly important to **rapidly** and **accurately diagnose** life-threatening **cardiac emergencies** such as **acute coronary syndromes (ACS)**<sup>(10)</sup> and **acute heart failure (HF)**<sup>(11)</sup>. However, diagnostic difficulties may lead to “over-admission” to in-hospital patient care for patients with suspected ACS with a negative effect on costs and resource utilization<sup>(12)</sup>.

**Biomarker tests** have been shown to contribute to **efficient triage** and improved patient management of acute cardiac conditions<sup>(13)</sup>. Cardiac biomarkers with a high negative predictive value (**NPV**) allow **rapid discharge** from the ED, whereas biomarkers with a high positive predictive value (**PPV**) are useful for **risk stratification** and **therapy guidance**.


This booklet describes the use of **cardiac markers** in the **diagnosis and risk stratification of ED patients** with signs and symptoms of cardiac emergencies such as ACS and acute HF. Emphasis is placed on **cardiac troponin** and **B-type natriuretic peptides** with reference to the most recent evidence-based professional recommendations<sup>(6,14)</sup>.

## OUR SPECIAL THANKS GO TO

Dr Pierre-Frédéric Keller

Cardiology Department / Intensive Care Unit  
University Hospitals - Geneva, Switzerland

**for his comprehensive review of this booklet.**

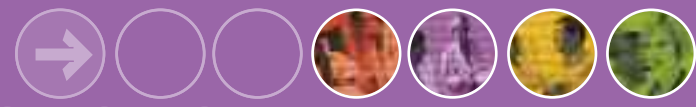

# Acute Coronary Syndromes



## Definition and classification

**Acute coronary syndromes (ACS)** refer to a constellation of clinical symptoms as a result of acute myocardial ischemia. It may range from a potentially reversible phase (**unstable angina**) to irreversible cell death (**myocardial infarction**).

The diagnosis and risk stratification of ACS is based on the integration of:

- the patient's presenting **symptoms**,
- **ECG** abnormalities,
- measurement of a **biomarker** of cardiac necrosis (**cardiac troponin**).

This results in the distinction of **three categories** (Table 1):

- unstable angina (**UA**),
- non-ST-segment elevation myocardial infarction (**NSTEMI**),
- ST-segment elevation myocardial infarction (**STEMI**).

The distinction between ACS categories is clinically important and drives the decision for type and intensity of therapeutic intervention.

- **ECG**: identifies approximately 1/3<sup>rd</sup> of ACS patients with persistent ST-segment elevation (**STEMI**) who require immediate reperfusion.
- **Cardiac troponin**: distinguishes the 2/3<sup>rd</sup> of ACS patients without ST-segment elevation (**NSTEACS**; non-ST-segment elevation ACS) who require either a conservative (**UA**) or early-invasive approach (**NSTEMI**).



**Table 1: Distinguishing features of acute coronary syndromes** <sup>(14, 15)</sup>

	MYOCARDIAL INFARCTION		
	NSTEMACS		STEMI
	Unstable Angina	NSTEMI	
<b>Pathophysiology</b>	Ischemia without necrosis	Ischemia with necrosis	
	Partially or transiently obstructive thrombus		Complete obstruction by intracoronary thrombus
<b>Clinical features</b> • Physical examination and history	Chest pain (angina and associated features) and presence of risk factors		
• Typical presenting symptoms	Severe angina (new onset, crescendo or rest angina)	Prolonged “crushing” chest pain, more severe and wider radiation than usual angina	
<b>12-lead ECG*</b>	No abnormalities, transient ST-elevation, ST-depression or T-wave inversion		Persistent ST-elevation, new left bundle branch block (LBBB)
<b>Cardiac troponin</b> Measurement on arrival and at 6 h	Negative (2x)	Positive	Positive**
<b>Therapeutic intervention</b>	Non-invasive (conservative)	Early-invasive	Immediate reperfusion

\* Observation of dynamic profiles is more informative (repeat or continuous monitoring).

\*\* Useful for confirmation, but availability of cTn test result should not delay therapeutic intervention.



## Pathophysiology

UA, NSTEMI and STEMI have a common pathophysiological origin related to **atherosclerotic coronary artery disease (CAD)** <sup>(15, 16)</sup>. Progression of the atherosclerotic plaque may result in its erosion or rupture with subsequent activation of blood platelets and coagulation factors leading to the formation of an intracoronary thrombus. Intracoronary obstruction results in loss of blood flow to the myocardium, causing **ischemia** (imbalance between oxygen supply and demand) and ultimately myocardial death (**necrosis**).



In rare cases, ACS may occur due to ischemia in the absence of occlusive atherosclerosis (e.g. coronary spasm such as Prinzmetal's angina, cocaine abuse, or coronary artery inflammation such as Kawasaki disease) <sup>(15, 16)</sup>.



## Signs and symptoms

The chief complaint in ACS is **angina**, defined as central **chest pain** or discomfort that occurs due to inadequate delivery of oxygen to the heart muscle. The pain may radiate to the neck, jaw or (left) arm. "Typical stable angina" is provoked by exertion or emotional stress and is relieved by rest or sublingual nitrates, but not in the case of unstable angina or myocardial infarction.

Often, the discomfort is diffuse (not localized, not positional) and may be accompanied by sweating, dyspnea, nausea, syncope. Chest pain is not specific for ACS and may also occur in other cardiac and non-cardiac conditions (Table 2). Among the cardiac conditions, aortic dissection and pericarditis must be excluded before initiation of therapy for myocardial infarction.

**Table 2:** Causes of chest pain or discomfort

CARDIAC	NON-CARDIAC
<ul style="list-style-type: none"><li>• <b>Acute coronary syndrome</b></li><li>• <b>Aortic dissection</b></li><li>• <b>Pericarditis</b></li><li>• <b>Myocarditis</b></li><li>• <b>Valvular disease</b></li></ul>	<ul style="list-style-type: none"><li>• <b>Gastrointestinal</b> Esophageal spasm or reflux Peptic ulcer</li><li>• <b>Pulmonary</b> Pneumonia Pulmonary embolism Pneumothorax</li><li>• <b>Neurological</b> (nerve root pain, herpes zoster)</li><li>• <b>Musculoskeletal</b> (e.g. osteochondritis)</li></ul>


**Atypical** presentations of ACS (e.g. epigastric pain, upper back pain or dizziness) are frequent in younger (25-40 years) and older (>75 years) patients and in women. Asymptomatic myocardial ischemia (**silent ischemia**) is particularly common among diabetics.



## Electrocardiogram

A number of typical ECG abnormalities are often observed during ACS, including **ST-segment depression** or **T-wave inversion** in UA/NSTEMI and **ST-segment elevation** in the early phase of STEMI (Figure 1).

The ECG may also provide information on the anatomic location, extent and severity of the coronary lesion as well as the presence of complications of acute MI with a poor prognosis (e.g. arrhythmias and conduction abnormalities such as bundle branch block and heart block)<sup>(17)</sup>.

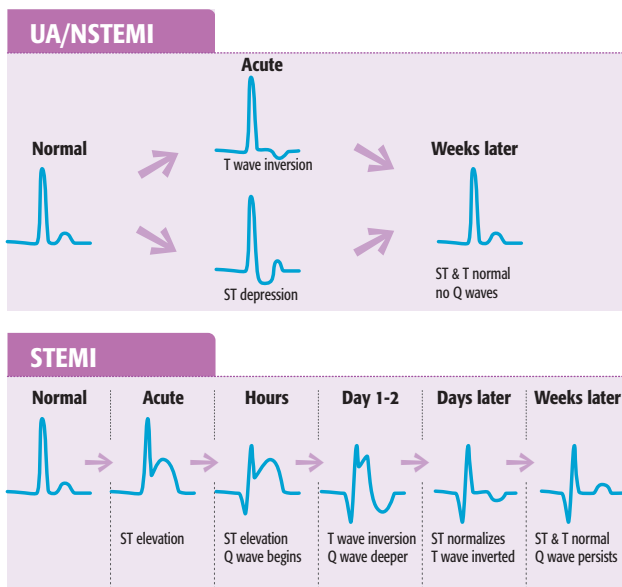

**Figure 1:** ECG evolution in ACS<sup>(18)</sup>

ST-segment elevation is the hallmark finding in STEMI, but may also be observed in other conditions (e.g. acute pericarditis and left ventricular hypertrophy). Abnormal Q-waves are less frequently observed nowadays due to earlier interventions.



## Diagnosis and risk stratification

### DECISION-MAKING ALGORITHM

The objective of the initial evaluation of patients with chest pain and suspected ACS is to address the following two questions<sup>(14, 15)</sup>:

#### 1. Differential diagnosis:

What is the likelihood that the patient's symptoms represent ACS due to underlying coronary artery disease?

#### 2. Risk stratification:

What is the likelihood that the patient will experience an adverse cardiovascular outcome (e.g. death, myocardial infarction, recurrent ischemia, stroke, heart failure)?

The outcome of this assessment (Figure 2) provides information that guides decision-making in terms of **selection, timing and intensity of therapeutic intervention, patient disposition** (coronary care unit, high dependency unit, etc.) and further investigations (**stress testing, angiography, etc.**).

Risk stratification can be refined by integration of cardiac necrosis markers and clinical factors into a **risk score**. Examples of validated risk scores are the **TIMI** and **GRACE** risk scores<sup>(19)</sup>.

Calculators for these scores are available on-line:

[www.timi.org](http://www.timi.org)


[www.outcomes.org/grace](http://www.outcomes.org/grace)







# Acute Coronary Syndromes

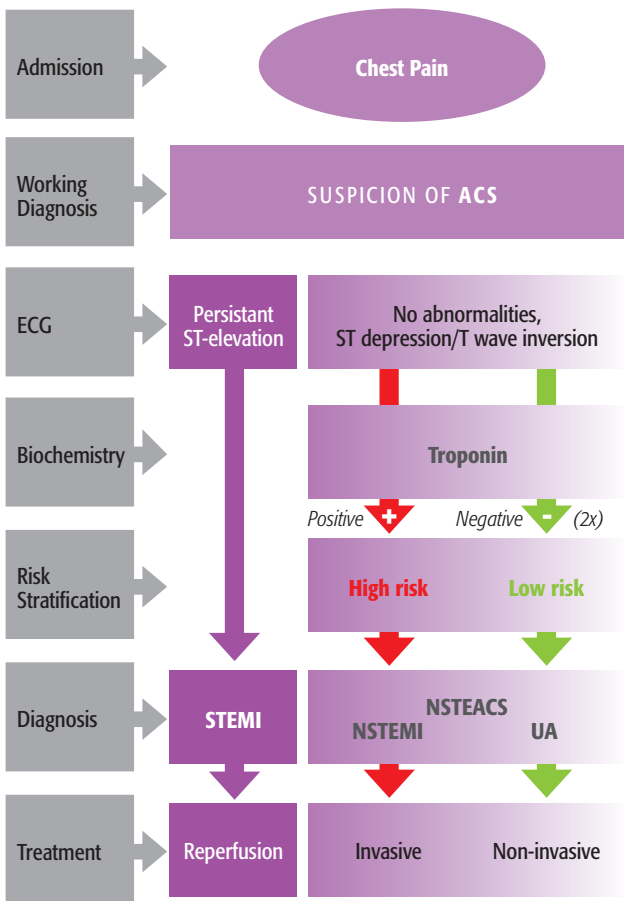

**Figure 2:** Approach to diagnosis and risk stratification (ESC Guidelines)<sup>(15)</sup>

**Serial measurement** of cTn (on arrival and after 6 hours) is required in patients without ST-elevation on ECG.



## Cardiac necrosis markers

**Myocardial infarction** is defined as myocardial cell death (**necrosis**) due to prolonged ischemia. Myocardial necrosis is recognized by the appearance in blood of different proteins that are released from the damaged myocytes. The best described and most widely available biomarkers of myocardial necrosis include cardiac troponin I and T (**cTnI, cTnT**), the MB fraction of creatine kinase (**CK-MB**) and **myoglobin**. These cardiac necrosis markers show important differences in key properties such as diagnostic performance (Table 3) and kinetic profile (Figure 3).


Because recognition of acute MI is important for prognosis and therapy selection, measurement of cardiac necrosis markers is indicated in all patients with suspected ACS<sup>(14, 15, 16)</sup>:

- **Cardiac troponin is the preferred cardiac necrosis biomarker.**
- CK-MB is an acceptable alternative when cTn is not available.

**Table 3:** Properties of cardiac necrosis markers<sup>(14)</sup>

MARKER	CARDIAC SPECIFICITY	TEMPORAL PROFILE		
		TIME TO FIRST DETECTION	MEAN TIME TO PEAK ELEVATION	DURATION OF ELEVATION
Myoglobin	+	1-3 h	6-7 h	12-24 h
CK-MB	+++	3-4 h	24 h	24-36 h
cTnI	++++	3-6 h	24 h	5-10 days
cTnT	++++	3-6 h	24 h	5-14 days

# Acute Coronary Syndromes

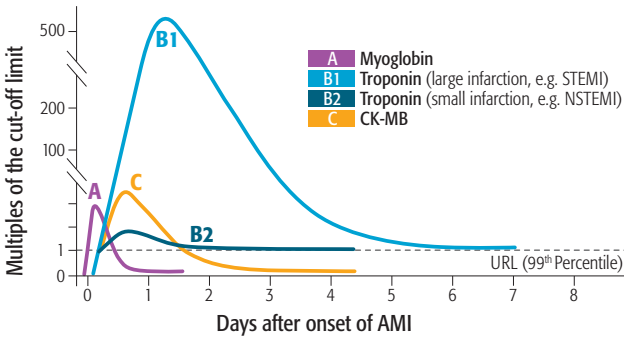

**Figure 3:** Temporal profile of cardiac necrosis markers after acute myocardial infarction <sup>(13)</sup>

Biomarker concentrations are plotted as multiples of the cut-off for AMI, i.e. any measurement exceeding the 99<sup>th</sup> percentile of a normal reference population (URL = upper reference limit). Troponin shows small elevations above the URL in small infarctions (typically in NSTEMI) but may rise to 20 to 50 times the URL with large infarctions (typically in STEMI).

## CLINICAL UTILITY

### ADVANTAGE

High sensitivity and NPV.  
**Early detection of MI** (early rule-out) and detection of **reperfusion**.  
 Detection of **reinfarction**.  
 Large clinical experience, previous "gold standard" for myocardial necrosis (best alternative if cTn assays are not available).  
 Superior **sensitivity** and **specificity**.  
**Current biomarker of choice for detection of myocardial injury**.  
 Powerful tool for risk stratification and therapy selection.  
 Detection of recent MI up to 2 weeks.

### DISADVANTAGE

Low specificity in presence of skeletal muscle injury and renal insufficiency.  
 Rapid clearance.  
 Reduced specificity in presence of skeletal muscle injury.  
 Gender-specific cut-off values.  
 Not an early marker of myocardial necrosis; serial testing needed when first result is normal.  
 Not an early marker of myocardial necrosis; serial testing needed when first result is normal.  
 Reduced ability to discriminate reinfarction (serial testing needed).

CK-MB: creatine kinase MB fraction; cTn: cardiac troponin




## Cardiac troponin

BIOMARKER OF CHOICE FOR DETECTING CARDIAC INJURY

The troponin complex consists of 3 subunits (I, T and C) and is essential for the regulation of skeletal and cardiac muscle contraction (Figure 4). In contrast to troponin C, **cardiac-specific isoforms of troponin T and I exist**. Specific antibodies have been raised against these cardiac isoforms, and form the basis of widely available quantitative and reliable cardiac troponin (cTn) assays.

Due to its **high sensitivity** and nearly absolute myocardial tissue **specificity**, **cTn assays** have become the **cornerstone** in the **diagnosis** and **risk stratification** of ACS <sup>(14, 15, 16)</sup>.

### CARDIAC MYOCYTE

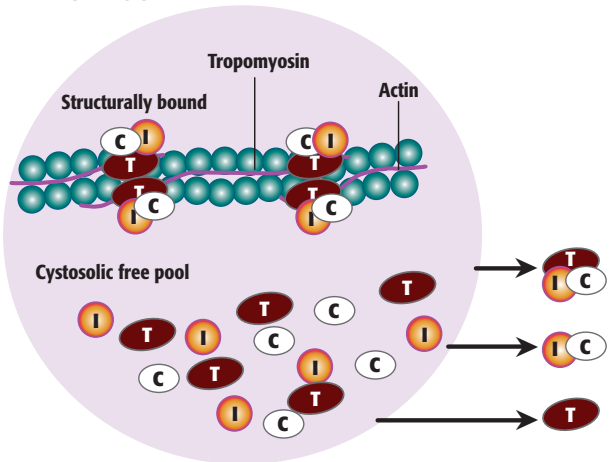

**Figure 4:** Troponin is a marker of cardiac damage <sup>(20)</sup>

The troponin complex is essential for the calcium-mediated regulation of muscle contraction. It consists of 3 subunits (troponin I, T and C) structurally bound on the actin filament. The cytosol of the cardiac myocyte contains unbound troponins which are released into the circulation upon damage.

**Elevated cTn is only diagnostic for acute MI when there is evidence of myocardial ischemia** (i.e. clinical symptoms and/or ECG abnormalities or imaging evidence) <sup>(16)</sup>.





## Acute Coronary Syndromes

The utility of cTn in the diagnosis of acute MI is contingent on the use of the proper **cut-off** value and the **timing** of measurement (Table 4).

**Table 4:** Diagnosis of acute MI: considerations for use of cTn<sup>(14, 16)</sup>

CRITERION	RECOMMENDATION
<b>Decision level for MI</b>	The <b>99<sup>th</sup> percentile*</b> of a normal reference population (URL=upper reference limit)**.
<b>Assay precision</b>	The optimal precision at the URL should be <b>≤ 10%</b> (total CV)*.
<b>Timing of measurement</b>	<b>Admission and 6 h later</b> At 12-24 h (if earlier measurements are normal and clinical suspicion of MI is high)

\* See “Frequently Asked Questions” section, page 22, for further information.

\*\* The actual value of the URL depends on the particular assay that is being used. The 1xURL cut-off applies to MI due to a primary coronary event, the most common form. Higher cut-off levels are indicated for MI associated with percutaneous coronary intervention (PCI; cut-off is 3xURL) and coronary artery bypass grafting (CABG; cut-off is 5xURL).

Elevated cTn is specific for cardiac damage, but not for coronary disease and can also be elevated **in the absence of ACS** (Table 5).

**Table 5:** Causes of elevated cTn in absence of ACS<sup>(16, 20)</sup>

<b>Demand ischemia</b>	<ul style="list-style-type: none"><li>• Critically ill patients (respiratory failure, sepsis)</li><li>• Arrhythmias</li><li>• Aortic dissection</li></ul>
<b>Myocardial ischemia</b>	<ul style="list-style-type: none"><li>• Hemorrhagic or ischemic stroke</li></ul>
<b>Direct myocardial damage</b>	<ul style="list-style-type: none"><li>• Trauma including cardiac surgery, ablation, pacing, etc.</li><li>• Rhabdomyolysis</li><li>• Drug toxicity (e.g. chemotherapy)</li><li>• Cardiac infiltrative diseases</li><li>• Inflammatory diseases (myocarditis, pericarditis)</li></ul>
<b>Myocardial strain</b>	<ul style="list-style-type: none"><li>• Heart failure</li><li>• Pulmonary embolism</li><li>• Severe pulmonary hypertension</li><li>• Extreme exertion</li></ul>
<b>Renal failure</b>	



# Acute Heart Failure



## Definition and classification

Heart failure (HF) is a complex clinical syndrome in which the pumping function of the heart becomes insufficient (**ventricular dysfunction**) to meet the needs of the vital systems and tissues of the body. A unified and practical definition of HF has recently been put forward by the ESC (Table 6).

Many **descriptive terms** are used to characterize and classify patients with HF:

- **TEMPORAL PRESENTATION**

**Acute** or stable **chronic** HF.

- **LEFT VENTRICULAR EJECTION FRACTION (LVEF)**

**Systolic** (LVEF < 40%) or **diastolic** HF (i.e. HF with normal ejection fraction, LVEF > 40-50%). These are not strict distinct entities and most HF patients have evidence of both. Diastolic HF is more common in female and the elderly.

- **LOCATION**

**Right** HF (congestion of systemic veins causing peripheral edema or hepatomegaly) or **left** HF (most common, congestion of pulmonary veins causing pulmonary edema).

The **severity** of HF is commonly described by the **NYHA functional classification system**, based on symptoms and exercise capacity (Table 7).





**Table 6:** ESC definition of heart failure <sup>(21)</sup>

<b>Heart failure is a clinical syndrome in which patients have the following features:</b>	
<b>Symptoms typical of heart failure</b>	Breathlessness at rest or on exercise Fatigue Tiredness Ankle swelling
<b>AND</b>	
<b>Signs typical of heart failure</b>	Tachycardia Tachypnea Pulmonary rales Pleural effusion Raised jugular venous pressure Peripheral edema Hepatomegaly
<b>AND</b>	
<b>Objective evidence of a structural or functional abnormality of the heart at rest</b>	Cardiomegaly Third heart sound Cardiac murmurs Abnormality on the echocardiogram Raised natriuretic peptide concentration

**Table 7:** Severity of heart failure: the NYHA classification <sup>(21)</sup>

<b>NYHA CLASS</b>	<b>DESCRIPTION</b>
<b>Class I</b>	No limitation of physical activity; ordinary physical activity does not cause undue fatigue, palpitation or dyspnea.
<b>Class II</b>	Slight limitation of physical activity; comfortable at rest, but ordinary physical activity results in fatigue, palpitation or dyspnea.
<b>Class III</b>	Marked limitation of physical activity; comfortable at rest, but less than ordinary physical activity results in fatigue, palpitation or dyspnea.
<b>Class IV</b>	Unable to carry on any physical activity without discomfort; symptoms at rest, if any physical activity is undertaken, discomfort is increased.

NYHA: New York Heart Association



## Pathophysiology

HF is a progressive and chronic disease, worsening over time. HF is caused by progressive **remodeling** of the heart, a process that changes its size and shape and subsequently **impairs the function of the ventricles**:

- **Systolic HF**: thinning and weakening of the ventricle walls results in dilation and a reduced capacity to eject blood (reduced ejection fraction).
- **Diastolic HF**: thickening and stiffening of the ventricles due to hypertrophy results in impaired relaxation (preserved ejection fraction).

Any structural or functional disorder that leads to deterioration of heart muscle function may lead to HF. Coronary heart disease is the initiating cause in about 70% of patients with HF <sup>(21)</sup>.

**Acute HF** is either due to **acute decompensation** of previously stable chronic HF (ADCHF, 63%) or new onset HF (37%) <sup>(22)</sup>. The **causes** and **precipitating factors** are listed in Table 8. ACS is the most frequent precipitating factor of new onset acute HF, whereas non-compliance with therapy is the major cause of ADCHF <sup>(22)</sup>.

**Table 8:** Causes and precipitating factors of acute heart failure <sup>(21)</sup>


<b>Ischemic heart disease</b>	ACS Complications of acute MI
<b>Valvular diseases</b>	Valve stenosis Regurgitation Endocarditis Aortic dissection
<b>Myopathies</b>	Post-partum cardiomyopathy Acute myocarditis
<b>Hypertension/arrhythmia</b>	
<b>Circulatory failure</b>	Sepsis Anemia Pulmonary embolism
<b>Decompensation of pre-existing chronic HF</b>	Lack of compliance with treatment Volume overload Infections (especially pneumonia) Cerebrovascular insult Surgery Renal dysfunction Asthma, COPD Drug abuse, alcohol abuse





Signs and symptoms

- Patients with acute HF have an **altered hemodynamic profile** manifested by signs of **congestion** (“wet”) and **hypoperfusion** (“cold”) (Figure 5).
- **Dyspnea** is the most common presenting symptom, but this is not unique for acute HF (Table 9).
- Frequent **co-morbidities** in acute HF include coronary heart disease, hypertension, atrial fibrillation, diabetes, valvular disease, anemia, COPD and renal failure<sup>(22)</sup>.

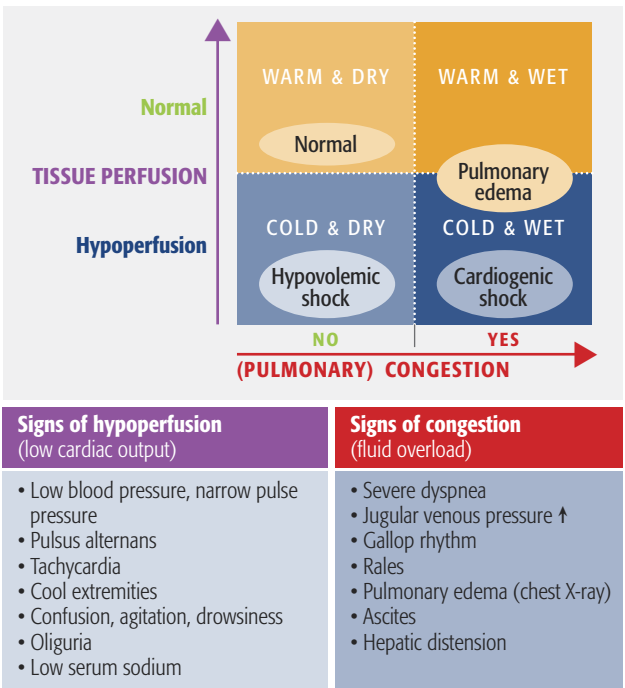

**Figure 5:** Hemodynamic profile in acute HF<sup>(21)</sup>

Clinical classification according to the hemodynamic profile (modified Forrester classification).



**Table 9:** Causes of dyspnea

CARDIOVASCULAR	NON-CARDIOVASCULAR
<b>CARDIAC</b> <ul style="list-style-type: none"> <li>• Heart failure</li> <li>• ACS</li> <li>• Significant valvular disease</li> <li>• Arrhythmias (especially atrial fibrillation)</li> <li>• Constrictive pericarditis/cardiac tamponade</li> <li>• Restrictive cardiomyopathy</li> </ul>	<b>RESPIRATORY</b> <ul style="list-style-type: none"> <li>• Pneumonia</li> <li>• Asthma</li> <li>• COPD</li> <li>• Pneumothorax</li> <li>• Pleural effusion</li> <li>• Upper airway obstruction</li> <li>• Pneumonitis/pulmonary fibrosis</li> </ul>
<b>NON-CARDIAC</b> <ul style="list-style-type: none"> <li>Pulmonary embolism</li> <li>Pulmonary hypertension</li> </ul>	<b>OTHERS</b> <ul style="list-style-type: none"> <li>• Anemia</li> <li>• Thyrotoxicosis</li> <li>• Metabolic, e.g. acidosis</li> <li>• Chest wall pain (pleuritic/musculoskeletal)</li> <li>• Skeletal abnormalities</li> <li>• Neuromuscular (diaphragmatic weakness)</li> <li>• Anxiety/psychogenic</li> </ul>



## Diagnosis and risk stratification

### DECISION-MAKING ALGORITHM

Diagnosing acute heart failure is difficult because there are many varied and often **non-specific clinical symptoms** (dyspnea, chest pain, fatigue, cough). The diagnosis, therefore, is based on the combined use of patient history and physical examination, ECG, chest X-ray, echocardiography and laboratory tests, including **natriuretic peptides** (Figure 6) <sup>(21)</sup>.

An **accurate** and **early diagnosis** is essential to target appropriate therapy and improve patient outcome. Apart from determining the presence and type of HF, the medical decision-making process (risk stratification and therapy selection) is further based on a complete understanding of underlying etiology, hemodynamic profile and stage/severity <sup>(21)</sup>.



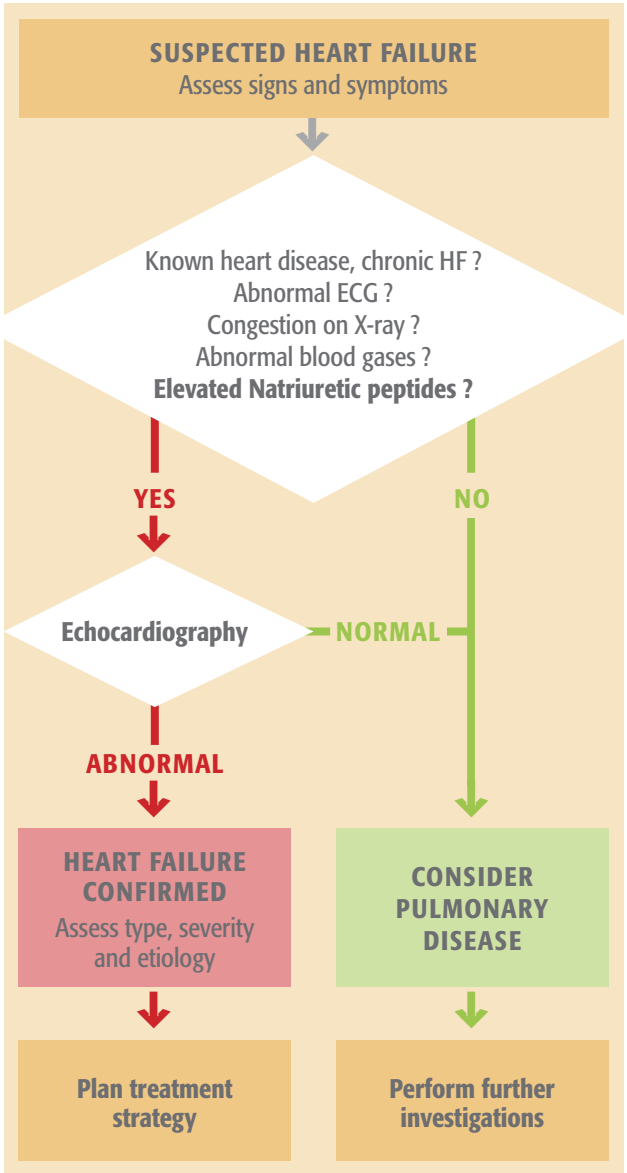


Figure 6: Evaluation of suspected acute HF (ESC Guidelines) <sup>(21)</sup>



## Natriuretic peptides

### BIOMARKERS OF MYOCYTE STRESS

The natriuretic peptide family consists of 3 peptides: atrial natriuretic peptide (ANP), brain (or B-type) natriuretic peptide (BNP) and C-type natriuretic peptide (CNP). These neurohormones are released in response to hemodynamic stress and are involved in the regulation of intravascular volume homeostasis<sup>(23, 24)</sup>.

**BNP** is secreted by the ventricles, and to a lesser extent by the atria, and appears in blood after cleavage of the precursor molecule proBNP. This cleavage also results in the release of **NT-proBNP**, the N-terminal counterpart (Figure 7). Therefore, the blood levels of both molecules are increased in HF.

BNP and NT-proBNP assays show **similar clinical performance characteristics** and their levels are well correlated<sup>(25)</sup>.


- Different cut-off values of BNP and NT-proBNP are in use.
- Absolute levels of BNP and NT-proBNP are not interchangeable.

**Table 10:** Properties of B-type natriuretic peptides<sup>(24)</sup>

CHARACTERISTIC	BNP	NT-proBNP
<b>Biologically active</b>	Yes	No
<b>Prohormone fragment</b>	C-terminal (proBNP 77-108) 32 amino acids	N-terminal (proBNP 1-76) 76 amino acids
<b>Half-life (min)</b>	20	60-120
<b>In vitro sample stability</b> (room temperature)	4 hours	> 3 days
<b>Sample type</b>	Whole blood, plasma (EDTA)	Plasma (heparin) or serum
<b>Assay measuring range (pg/mL)</b>	5 – 5,000	20 – 25,000*

\* VIDAS® NT-proBNP (bioMérieux)



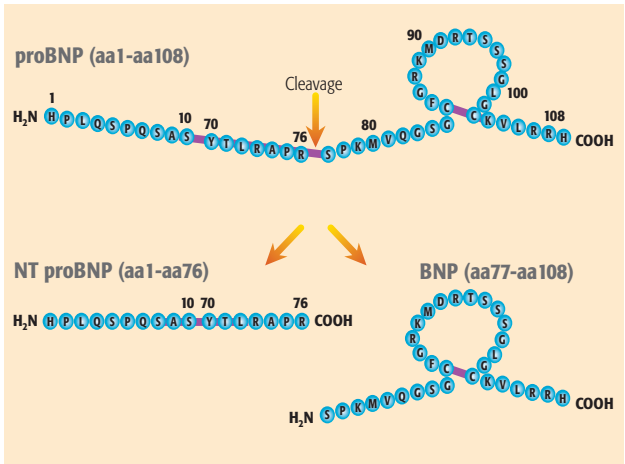

**Figure 7:** Release of BNP and NT-proBNP <sup>(24)</sup>

BNP and NT-proBNP are quantitative markers of cardiac stress that are released into blood after cleavage of the precursor protein proBNP.




## → NT-proBNP for evaluation of HF in patients with acute dyspnea

The measurement of B-type natriuretic peptides is recommended in many guidelines as an integral part in the HF diagnostic work-up<sup>(6,21,25)</sup>. This is particularly useful in the assessment of patients with acute dyspnea in the **ED setting**, where NT-proBNP has utility as both a **rule-out** and a **rule-in** test for HF (Figure 8).

**NT-proBNP** is highly **sensitive** and **specific** for the **diagnosis** or **exclusion** of **acute HF** and is a powerful and cost-effective adjunctive tool for the clinician in the diagnosis and triage of patients with acute dyspnea<sup>(26)</sup>.

The NT-proBNP level must be interpreted in conjunction with a thorough history and physical examination of the patient, as it can also be elevated in the absence of HF (Table 11).

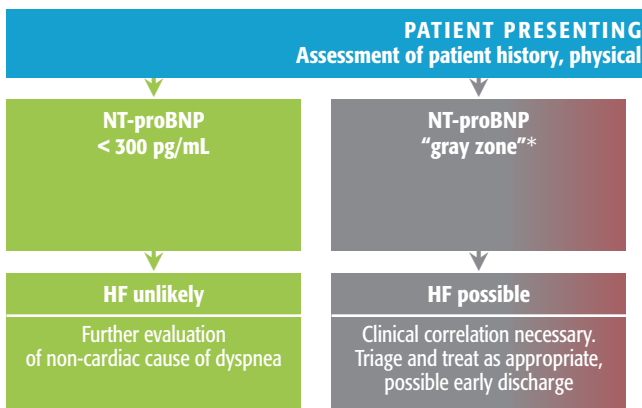


**Figure 8:** NT-proBNP in the evaluation and triage of ED patients with acute dyspnea<sup>(26)</sup>

In this cost-effective and clinically validated algorithm (International NT-proBNP Consensus Panel), a **single cut-off value of 300 pg/mL** is used for **exclusion**, whereas age-adjusted cut-off values are used to rule-in HF.

**Table 11:** Elevated NT-proBNP in absence of HF<sup>(27)</sup>

- **Heart muscle disease**  
(e.g. myocarditis, cardiomyopathies, amyloidosis)
- **Heart valvular disease**  
(e.g. aortic/mitral stenosis and regurgitation)
- **Arrhythmia** (atrial fibrillation)
- **Acute coronary syndrome**
- **Stroke**
- **Pulmonary embolism**
- **Chronic pulmonary disease**  
(e.g. COPD, pulmonary artery hypertension)
- **Anemia**
- **Renal failure**
- **Diabetes mellitus**
- **Critical illness** (e.g. sepsis, burns, ARDS)

**WITH ACUTE DYSPNEA**  
examination, ECG, chest X-ray and NT-proBNP

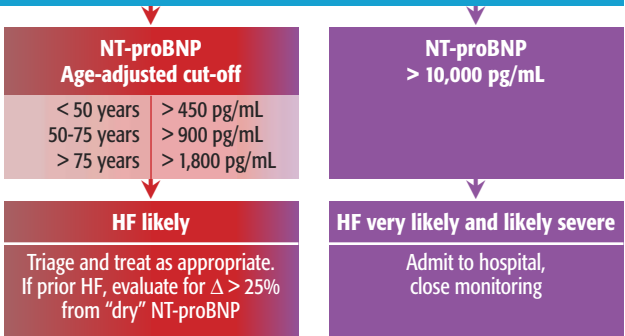

\* The area between the rule-out (< 300 pg/mL) and the rule-in (age-adjusted) cut-off values is designated as the "gray zone".

A "gray zone" NT-proBNP value **should not be ignored** because it is associated with **worse outcomes** compared with NT-proBNP values below the rule-out cut-off<sup>(28)</sup>. See Frequently Asked Questions section, page 25, for further information.

# Frequently Asked Questions


## Cardiac troponin and ACS

### What is the optimal cTn decision limit for myocardial infarction (i.e. what is the difference between the 99<sup>th</sup> percentile of a normal reference population and the 10% CV limit) ?

According to the universal definition of myocardial infarction, the 99<sup>th</sup> percentile of a normal reference population is recommended as the **decision limit** <sup>(16)</sup>. Furthermore, the optimal precision (total coefficient of variation) at this level should be  $\leq 10\%$  <sup>(16)</sup>.

For most cTn assays, however, the lowest concentration with a total CV of 10% is above the 99<sup>th</sup> percentile reference. Therefore, experts have in the past proposed to use the higher 10% CV value as an acceptable alternative decision limit <sup>(29)</sup>. This view is being abandoned and, irrespective of the total CV at the 99<sup>th</sup> percentile value, only the 99<sup>th</sup> percentile value is recommended as the MI decision limit <sup>(30)</sup>.

**Intermediate cTn values** (i.e. above the 99<sup>th</sup> percentile but below the 10% CV limit) have **prognostic value** and it is not appropriate to disregard them (Figure 9) <sup>(31, 32)</sup>.

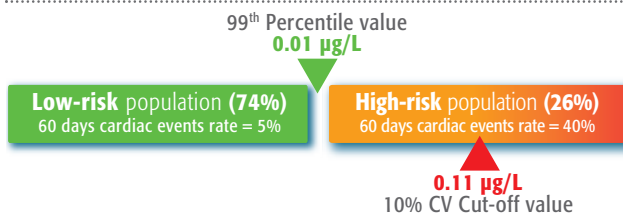

**Figure 9:** VIDAS<sup>®</sup> Troponin I Ultra (bioMérieux) 99<sup>th</sup> percentile and prediction of cardiac events <sup>(32)</sup>

In a cohort of suspected ACS patients (n=302), cTnI values were measured on admission on the VIDAS system. Patients with cTnI above the 99<sup>th</sup> percentile (26%) had a 9-fold increased risk for adverse cardiac events (MI or death) at 60 days <sup>(32)</sup>.





## What is the clinical impact of the new generation of ultra-sensitive cTn assays?

The new generation of cTn assays allows the detection of **lower concentrations with higher precision** <sup>(33)</sup>. With the newer assays, cTn elevations are detected **within about 2 hours from the onset of symptoms**. The **earlier detection** of myocardial necrosis by the newer cTn assays could result in **decreased costs** and **better patient care** because of earlier triage to an invasive strategy or earlier discharge <sup>(34)</sup>.

On the other hand, assays with improved sensitivity will also lead to an increase in the number of patients with slightly elevated cTn concentrations due to “minor” myocardial damage, not always related to an ischemic cause. Yet, there is increasing evidence that these minor previously undetectable elevations of cTn are prognostic in a variety of populations, including the general community, and patients with heart failure and/or stable coronary artery disease <sup>(34)</sup>.

## Can a false-negative cTn result be expected?

To diagnose myocardial infarction, it is important to perform **serial measurement** of cTn because the first measurement on admission may still be negative when the **patient presents early after the onset of symptoms** <sup>(14, 16)</sup>. Analytical false-negatives may occur in the presence of **circulating troponin autoantibodies** <sup>(35)</sup>.

## Can a false-positive cTn result be expected?

Analytical false-positives may occur due to the presence of potentially interfering factors such as **heterophilic antibodies** or **rheumatoid factor** <sup>(20)</sup>. Medical false-positives in the absence of ACS may occur with many **underlying diseases** (Table 5). However, the term “false-positives” is misleading because, **even without cardiac ischemia, elevated cTn correlates with a worse prognosis** <sup>(36)</sup>.

Serial cTn testing is necessary to distinguish an acute coronary event, characterized by a typical rising and falling pattern, from other more chronic conditions <sup>(33)</sup>. Changes of  $\geq 20\%$  in cTn are considered significant <sup>(14, 16)</sup>.



## What is the difference between cTnI and cTnT?

Both cTnI and cTnT are **equally effective** in the **diagnosis and prognosis of ACS**<sup>(14, 16)</sup>. There is only one manufacturer of cTnT assays, whereas cTnI assays are available from many sources with differences in standards and antibodies. Consequently, comparability of cTnI test results between methods is not optimal<sup>(37)</sup>. There is a different clearance pattern and the duration of cTnT elevation is longer (2 weeks) than cTnI (1 week)<sup>(14)</sup>. Renal failure leads more frequently to elevations in cTnT than cTnI<sup>(38)</sup>.

## Is cTn measurement useful in ICU patients?

Cardiac conditions are common in the ICU and cardiac biomarkers may play a role in providing **additional information in critically ill patients**<sup>(39)</sup>. Elevated cTn is frequently observed in such patients and is associated with **increased mortality and ICU length of stay**<sup>(40)</sup>. Elevated cTn alone, however, is not diagnostic of MI in critically ill patients, because it is also elevated in conditions such as renal failure, sepsis, trauma, heart failure and inflammatory disorders (see Table 5, page 11).

## What is the value of combined measurement of cTn and NT-proBNP in ACS?

Risk stratification is less accurate in patients with suspected ACS and normal cTn levels, because it relies only on clinical and ECG variables. Because NT-proBNP provides **independent prognostic information** irrespective of cTn status, its measurement has been recommended in ACS patients<sup>(14, 41)</sup>. This is particularly important in patients with **suspected ACS and normal cTn**, where **NT-proBNP > 474 pg/mL** was able to discriminate individuals **at higher risk**<sup>(42)</sup>.



## NT-proBNP and HF

### What is the difference between BNP and NT-proBNP?

BNP and NT-proBNP refer to the C-terminal and N-terminal cleavage products of proBNP, respectively. Due to differences in clearance, NT-proBNP has a longer half-life<sup>(23)</sup>.





Although numerically different values are reported, BNP and NT-proBNP generally have **equivalent diagnostic and prognostic properties**<sup>(25)</sup>. In direct comparison studies, however, NT-proBNP was reported to have **better accuracy for mild HF**<sup>(43)</sup> and HF with preserved ejection fraction<sup>(44)</sup>.

The main difference resides in their (pre)analytical properties which may favor **NT-proBNP as a more convenient molecule** to work with in clinical laboratories<sup>(45)</sup>:

- **Good harmony of results**, because the various NT-proBNP assays on the market are based on the same antibodies.
- **Better sample stability** at different temperatures and wider dynamic range, owing to its longer half-life.
- **Excellent precision** on automated instruments.

### How important is NT-proBNP in the “gray zone”?

NT-proBNP values **between the rule-out and rule-in cut-off values** for HF in the ED are referred to as **intermediate or gray zone** values (see Figure 8). In the **ICON study**, the gray zone was observed in **17% of dyspneic patients**, 54% of whom had HF as the ultimate diagnosis<sup>(28)</sup>. Gray zone patients usually have **mild HF with fairly good short-term outcome**. Nevertheless, a gray zone NT-proBNP value should not be ignored because it is associated with worse outcomes compared with NT-proBNP values below the rule-out cut-off<sup>(28)</sup>.

### Are there confounders that may affect the usefulness of NT-proBNP for HF diagnosis in ED patients?

Co-morbid conditions such as **COPD, renal disease, obesity, diabetes mellitus** and **atrial fibrillation** are frequently observed in ED patients with dyspnea. In the absence of HF, these conditions may result in elevated or decreased (in case of obesity) NT-proBNP levels. The influence of these conditions on the diagnostic accuracy of NT-proBNP for HF has been investigated in the PRIDE and ICON studies<sup>(46, 47, 48, 49, 50)</sup>.

With the exception of atrial fibrillation, these conditions do not influence the diagnostic accuracy of NT-proBNP and no changes are required to the recommended age-adjusted cut-off values.



## Can false-negative or false-positive NT-proBNP results be expected?

Despite severe symptoms, NT-proBNP may not be **very high** when HF is due to a cause upstream from the left ventricle such as in **mitral stenosis or acute mitral regurgitation** <sup>(25)</sup>. NT-proBNP also remains relatively low in patients presenting with HF symptoms that develop abruptly (< 1 hour), a rare condition referred to as **“flash” pulmonary edema** <sup>(25)</sup>.

Although NT-proBNP is elevated in many other conditions (see Table 11, page 21), the term “false-positives” is misleading because, **even without HF, elevated NT-proBNP is linked to an adverse outcome** <sup>(27)</sup>.

## What is the optimal timing of NT-proBNP measurement?

Timing of the initial NT-proBNP measurement in ED patients with suspected HF is critical because delayed measurement has been reported to be associated with delays in treatment and an increase in hospital mortality <sup>(51)</sup>. Subsequent **serial measurement at 4, 12 and 24 hours from admission** is useful in confirming the diagnosis of acute HF <sup>(52)</sup>. Variations in NT-proBNP levels are also predictive of outcome after hospital discharge <sup>(53)</sup>.

## Can NT-proBNP be used to guide therapy?

Natriuretic peptide levels are **commonly reduced by treatment** with diuretics, ACE inhibitors, angiotensin receptor blockers, aldosterone antagonists and cardiac resynchronization therapy <sup>(25)</sup>. Also, changes of NT-proBNP over months are predictive of outcome of HF <sup>(54)</sup>. This suggests that NT-proBNP could be **useful to guide therapy in selected cases**. However, there is no consensus among experts (mixed results of small controlled trials) and tailoring therapy to achieve a target NT-proBNP level is not warranted at this time <sup>(6, 25)</sup>.

## Is there a value for NT-proBNP measurement in the ICU?

Since increased NT-proBNP is not specific for HF and may be influenced by a variety of cardiac and non-cardiac conditions commonly seen in the ICU, the diagnostic accuracy for left ventricular failure is reduced <sup>(55)</sup>. In the ICU, NT-proBNP may be useful as a **general marker for cardiac dysfunction**, in distinguishing cardiogenic and non-cardiogenic pulmonary edema, and as an **aid in the timing of extubation** <sup>(25, 56)</sup>.





### Is there a value for NT-proBNP measurement in primary care?

Detection of patients with **early stage HF** may allow the **initiation of preventive therapies**. However, these patients are difficult to diagnose because of the very insidious signs and symptoms. Although NT-proBNP holds promise for the detection of these largely asymptomatic patients, it is not recommended yet for wide-scale population screening<sup>(57)</sup>. Nevertheless, **screening of targeted symptomatic patient populations** (e.g. diabetes, high blood pressure, ...) is valuable.

Compared with NT-proBNP values in ED patients with acute HF, lower values are expected in chronic HF patients in the community.

Therefore, the following cut-off values are recommended in primary care<sup>(58)</sup>:

AGE (years)	NT-proBNP (pg/mL)	INTERPRETATION
< 75 ≥ 75	< 125 < 450	<b>HF unlikely</b> Further investigation of non-cardiac causes.
< 75 ≥ 75	≥ 125 ≥ 450	<b>Left ventricular dysfunction possible</b> Further investigations needed.

Consequently, in primary care patients, NT-proBNP is only a **rule-out test to exclude significant cardiac disease**.

### Is there a value for combined measurement of cTn and NT-proBNP in HF?

Both NT-proBNP and cTn can identify patients with HF at increased risk for an adverse outcome<sup>(6)</sup>. Elevated cTn is not uncommon in HF and provides **prognostic information in patients hospitalized with acute decompensated heart failure**<sup>(59)</sup>. The combination of NT-proBNP and cTn is even more powerful, and, in ED patients with acute HF, the highest mortality rate was observed when both markers were elevated<sup>(60)</sup>. The therapeutic approach to such patients has not been established at this time.



# References

1. Lopez AD, Mathers CD, Ezzati M, et al. Global and regional burden of disease and risk factors, 2001: Systematic analysis of population health data. *Lancet*. 2006; 367: 1747-57.
2. Rosamond W, Flegal K, Furie K, et al. Heart disease and stroke statistics – 2008 update. A report from the American Heart Association Statistics Committee and Stroke Statistics Subcommittee. *Circulation*. 2008; 117: e25-146.
3. Scholte op Reimer WJM, Gitt AK, Boersma E, Simoons ML (eds.). *Cardiovascular Diseases in Europe. Euro Heart Survey 2006*. European Society of Cardiology, Sophia Antipolis, France.
4. Leal J, Luengo-Fernández R, Gray A, et al. Economic burden of cardiovascular diseases in the enlarged European Union. *Eur Heart J*. 2006; 27: 1610-9.
5. Bleumink GS, Knetsch AM, Sturkenboom et al. Quantifying the heart failure epidemic: Prevalence, incidence rate, lifetime risk and prognosis of heart failure The Rotterdam Study. *Eur Heart J*. 2004; 25: 1614-9.
6. Tang WH, Francis GS, Morrow DA, et al. National Academy of Clinical Biochemistry Laboratory Medicine Practice Guidelines: Clinical utilization of cardiac biomarker testing in heart failure. *Circulation*. 2007; 116: e99-109.
7. Ekelund U, Nilsson HJ, Frigyesi A, Torffvit O. Patients with suspected acute coronary syndrome in a university hospital emergency department: An observational study. *BMC Emerg Med*. 2002; 2: 1.
8. Januzzi JL, Camargo CA, Anwaruddin S, et al. The N-terminal pro-BNP investigation of dyspnea in the emergency department (PRIDE) study. *Am J Cardiol*. 2005; 95: 948-54.
9. Nawar EW, Niska RW, Xu J. National Hospital Ambulatory Medical Care Survey: 2005 Emergency Department Summary. *Advance Data from Vital and Health Statistics*, Number 386, June 29, 2007. National Center for Health Statistics ([www.cdc.gov/nchs](http://www.cdc.gov/nchs)).
10. Pope JH, Aufderheide TP, Ruthazer R, et al. Missed diagnoses of acute cardiac ischemia in the emergency department. *N Engl J Med*. 2000; 342: 1163-70.
11. Peacock WF. Using the emergency department clinical decision unit for acute decompensated heart failure. *Cardiol Clin*. 2005; 23: 569-88.
12. Forberg JL, Henriksen LS, Edenbrandt L, Ekelund U. Direct hospital costs of chest pain patients attending the emergency department: a retrospective study. *BMC Emerg Med*. 2006; 6: 6.
13. Jaffe AS, Babuin L, Apple FS. Biomarkers in acute cardiac disease. *J Am Coll Cardiol*. 2006; 48: 1-11.
14. Morrow DA, Cannon CP, Jesse JL, et al. National Academy of Clinical Biochemistry Laboratory Medicine Practice Guidelines: Clinical characteristics and utilization of biochemical markers in acute coronary syndromes. *Circulation*. 2007; 115: e356-75.
15. Bassand JP, Hamm CW, Ardissino D, et al. Task Force for Diagnosis and Treatment of Non-ST-Segment Elevation Acute Coronary Syndromes of the European Society of Cardiology. Guidelines for the diagnosis and treatment of non-ST-segment elevation acute coronary syndromes. *Eur Heart J*. 2007; 28: 1598-660.
16. Thygesen K, Alpert JS, White HD et al. Joint ESC/ACC/AHA/WHF Task Force for the Redefinition of Myocardial Infarction. Universal definition of myocardial infarction. *Circulation*. 2007; 116: 2634-53.
17. Zimetbaum PJ, Josephson ME. Use of the electrocardiogram in acute myocardial infarction. *N Engl J Med*. 2003; 348: 933-40.
18. Gupta A, Sabatine MS, O’Gara PT, Lilly LS. Acute coronary syndromes. In: *Pathophysiology of heart disease. A collaborative project of medical students and faculty*, 3<sup>rd</sup> edition. Lilly LS (ed). Lippincott Williams & Wilkins. Philadelphia, 2003.
19. Yan AT, Yan RT, Tan M, et al. Risk scores for risk stratification in acute coronary syndromes: useful but simpler is not necessarily better. *Eur Heart J*. 2007; 28: 1072-8.
20. Korff S, Katus HA, Giannitsis E. Differential diagnosis of elevated troponins. *Heart*. 2006; 92: 987-93.





21. Dickstein K, Cohen-Solal A, Filippatos G, et al. Task Force for Diagnosis and Treatment of Acute and Chronic Heart Failure 2008 of European Society of Cardiology. ESC Guidelines for the diagnosis and treatment of acute and chronic heart failure 2008: the Task Force for the Diagnosis and Treatment of Acute and Chronic Heart Failure 2008 of the European Society of Cardiology. Developed in collaboration with the Heart Failure Association of the ESC (HFA) and endorsed by the European Society of Intensive Care Medicine (ESICM). *Eur Heart J.* 2008; 29: 2388-442.
22. Nieminen MS, Brutsaert D, Dickstein K, et al. EuroHeart Failure Survey II (EHFS II): A survey on hospitalized acute heart failure patients: description of population. *Eur Heart J.* 2006; 27: 2725-36.
23. Martinez-Rumayor A, Richards M, Burnett JC, Januzzi JL. Biology of the natriuretic peptides. *Am J Cardiol.* 2008; 101 (Suppl.): 3A-8A.
24. Daniels LB, Maisel AS. Natriuretic peptides. *J Am Coll Cardiol.* 2007; 50: 2357-68.
25. Maisel A, Mueller C, Adams K Jr, et al. State of the art: Using natriuretic peptide levels in clinical practice. *Eur J Heart Fail.* 2008; 10: 824-39.
26. Januzzi JL, Chen-Tournoux AA, Moe G. Amino-terminal pro-B-type natriuretic peptide testing for the diagnosis or exclusion of heart failure in patients with acute symptoms. *Am J Cardiol.* 2008; 101 (Suppl.): 29A-38A.
27. Baggish AL, van Kimmenade RR, Januzzi JL. The differential diagnosis of an elevated amino-terminal pro-B-type natriuretic peptide level. *Am J Cardiol.* 2008; 101 (Suppl.): 43A-48A.
28. van Kimmenade RR, Pinto YM, Bayes-Genis A, et al. Usefulness of intermediate amino-terminal pro-brain natriuretic peptide concentrations for diagnosis and prognosis of acute heart failure. *Am J Cardiol.* 2006; 98: 386-90.
29. Apple FS, Wu AH, Jaffe AS. European Society of Cardiology and American College of Cardiology guidelines for redefinition of myocardial infarction: How to use existing assays clinically and for clinical trials. *Am Heart J.* 2002; 144: 981-6.
30. Apple FS, Parvin CA, Buechler KF, et al. Validation of the 99th percentile cutoff independent of assay imprecision (CV) for cardiac troponin monitoring for ruling out myocardial infarction. *Clin Chem.* 2005; 51: 2198-200.
31. Kontos MC, Shah R, Fritz LM, et al. Implication of different cardiac troponin I levels for clinical outcomes and prognosis of acute chest pain patients. *J Am Coll Cardiol.* 2004; 43: 958-65.
32. Apple FS, Smith SW, Pearce LA, et al. Use of the bioMérieux VIDAS® Troponin I Ultra assay for the diagnosis of myocardial infarction and detection of adverse events in patients presenting with symptoms suggestive of acute coronary syndromes. *Clin Chim Acta.* 2008; 390: 72-5.
33. Wu AH, Jaffe AS. The clinical need for high-sensitivity cardiac troponin assays for acute coronary syndromes and the role for serial testing. *Am Heart J.* 2008; 155: 208-14.
34. White HD. Will new higher-precision troponins lead to clarity or confusion? *Curr Opin Cardiol.* 2008; 23: 202-5.
35. Eriksson S, Hellman J, Pettersson K. Autoantibodies against cardiac troponins. *N Engl J Med.* 2005; 352: 98-100.
36. Alcalai R, Planer D, Culhaoglu A, et al. Acute coronary syndrome vs. nonspecific troponin elevation: Clinical predictors and survival analysis. *Arch Intern Med.* 2007; 167:276-81.
37. Panteghini M. Standardization of cardiac troponin I measurements: the way forward? *Clin Chem.* 2005; 51: 1594-7.
38. Freda BJ, Tang WH, Van Lente F, et al. Cardiac troponins in renal insufficiency: Review and clinical implications. *J Am Coll Cardiol.* 2002; 40:2065-71.
39. McLean AS, Huang SJ, Salter M. Bench-to-bedside review: The value of cardiac biomarkers in the intensive care patient. *Crit Care.* 2008; 12: 215.
40. Lim W, Qushmaq I, Devereaux PJ, et al. Elevated cardiac troponin measurements in critically ill patients. *Arch Intern Med.* 2006; 166: 2446-54.
41. Omland T, de Lemos JA. Amino-terminal pro-B-type natriuretic peptides in stable and unstable ischemic heart disease. *Am J Cardiol.* 2008; 101 (Suppl.): 61A-66A.
42. Weber M, Bazzino O, Navarro Estrada JL, et al. N-terminal B-type natriuretic peptide assessment provides incremental prognostic information in patients with acute coronary syndromes and normal troponin T values upon admission. *J Am Coll Cardiol.* 2008; 51:1188-95.
43. Emdin M, Passino C, Prontera C, et al. Comparison of brain natriuretic peptide (BNP) and amino-terminal proBNP for early diagnosis of heart failure. *Clin Chem.* 2007; 53: 1289-97.



44. O'Donoghue M, Chen A, Baggish AL, et al. The effects of ejection fraction on N-terminal ProBNP and BNP levels in patients with acute CHF: Analysis from the ProBNP Investigation of Dyspnea in the Emergency Department (PRIDE) study. *J Card Fail.* 2005; 11 (Suppl.): S9-14.
45. Ordóñez-Llanos J, Collinson PO, Christenson RH. Amino-terminal pro-B-type natriuretic peptide: Analytic considerations. *Am J Cardiol.* 2008; 101 (Suppl.): 9A-15A.
46. Tung RH, Camargo CA, Krauser D, et al. Amino-terminal pro-brain natriuretic peptide for the diagnosis of acute heart failure in patients with previous obstructive airway disease. *Ann Emerg Med.* 2006; 48: 66-74.
47. Anwaruddin S, Lloyd-Jones DM, Baggish A, et al. Renal function, congestive heart failure, and amino-terminal pro-brain natriuretic peptide measurement. Results from the ProBNP Investigation of Dyspnea in the Emergency Department (PRIDE) study. *J Am Coll Cardiol.* 2006; 47: 91-7.
48. Bayes-Genis A, Lloyd-Jones D, van Kimmenade RR, et al. Effect of body mass index on diagnostic and prognostic usefulness of amino-terminal pro-brain natriuretic peptide in patients with acute dyspnea. *Arch Intern Med.* 2007; 167: 400-7.
49. O'Donoghue M, Kenney P, Oestreicher E, et al. Usefulness of aminoterminal pro-brain natriuretic peptide testing for the diagnostic and prognostic evaluation of dyspneic patients with diabetes mellitus seen in the emergency department (from the PRIDE study). *Am J Cardiol.* 2007; 100: 1336-40.
50. Morello A, Lloyd-Jones DM, Chae CU, et al. Association of atrial fibrillation and amino-terminal pro-brain natriuretic peptide concentrations in dyspneic subjects with and without acute heart failure: Results from the ProBNP Investigation of Dyspnea in the Emergency Department (PRIDE) study. *Am Heart J.* 2007; 153: 90-7.
51. Maisel AS, Peacock WF, McMullin N, et al. Timing of immunoreactive B-type natriuretic peptide levels and treatment delay in acute decompensated heart failure: an ADHERE (Acute Decompensated Heart Failure National Registry) analysis. *J Am Coll Cardiol.* 2008; 52: 534-40.
52. Disomma S, Magrini L, Pittoni V, et al. Usefulness of serial assessment of natriuretic peptides in the emergency department for patients with acute decompensated heart failure. *Congest Heart Fail.* 2008; 14 (Suppl.): 21-4.
53. Bettencourt P, Azevedo A, Pimenta J, et al. N-terminal-pro-brain natriuretic peptide predicts outcome after hospital discharge in heart failure patients. *Circulation.* 2004; 110: 2168-74.
54. Masson S, Latini R, Anand IS, et al. Prognostic value of changes in N-terminal pro-brain natriuretic peptide in Val-HeFT (Valsartan Heart Failure Trial). *J Am Coll Cardiol.* 2008; 52: 997-1003.
55. Omland T. Advances in congestive heart failure management in the intensive care unit: B-type natriuretic peptides in evaluation of acute heart failure. *Crit Care Med.* 2008; 36 (Suppl.): S17-27.
56. Mueller C. The use of B-type natriuretic peptides in the intensive care unit. *Crit Care Med.* 2007; 35: 2438-9.
57. de Lemos JA, Hildebrandt P. Amino-terminal pro-B-type natriuretic peptides: testing in general populations. *Am J Cardiol.* 2008; 101 (Suppl.): 16A-20A.
58. Hildebrandt P, Collinson PO. Amino-terminal pro-B-type natriuretic peptide testing to assist the diagnostic evaluation of heart failure in symptomatic primary care patients. *Am J Cardiol.* 2008; 101 (Suppl.): 25A-28A.
59. Peacock WF IV, De Marco T, Fonarow GC, et al. Cardiac troponin and outcome in acute heart failure. *N Engl J Med.* 2008; 358: 2117-26.
60. Sakhuja R, Green S, Oestreicher EM, et al. Amino-terminal pro-brain natriuretic peptide, brain natriuretic peptide, and troponin T for prediction of mortality in acute heart failure. *Clin Chem.* 2007; 53:412-20.









# List of abbreviations

ACE	Angiotensin-converting enzyme
ACS	Acute coronary syndromes
ADCHF	Acute decompensation of chronic HF
ANP	Atrial natriuretic peptide
ARDS	Acute respiratory distress syndrome
BNP	Brain (or B-type) natriuretic peptide
CABG	Coronary artery bypass grafting
CAD	Coronary artery disease
CK-MB	Creatine kinase MB fraction
CNP	C-type natriuretic peptide
COPD	Chronic obstructive pulmonary disease
cTn	Cardiac troponin
CV	Coefficient of variation
ECG	Electrocardiogram
ED	Emergency department
ESC	European Society of Cardiology
GRACE	Global Registry of Acute Coronary Events
HF	Heart failure
ICU	Intensive care unit
LVEF	Left ventricular ejection fraction
MI	Myocardial infarction
NPV	Negative predictive value
NSTEACS	Non-ST-segment elevation acute coronary syndromes
NSTEMI	Non-ST-segment elevation myocardial infarction
NT-proBNP	N-terminal pro B-type natriuretic peptide
NYHA	New York Heart Association
PPV	Positive predictive value
SOB	Shortness of breath
STEMI	ST-segment elevation myocardial infarction
TIMI	Thrombolysis in myocardial infarction
UA	Unstable angina
URL	Upper reference limit

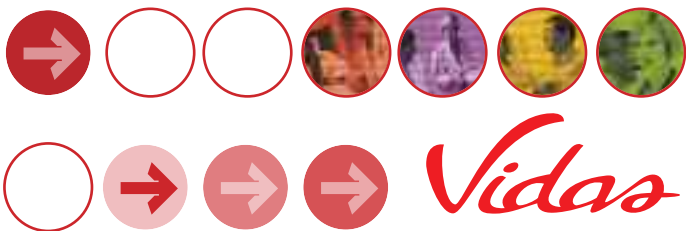

## Contribution of a Rapid Diagnostic Test Panel Patients presenting

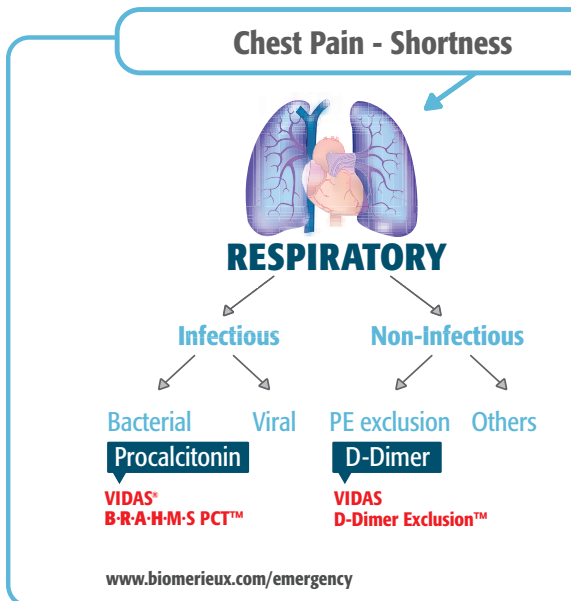

\* Four of the most frequent symptoms in patients presenting to Hospital ED.  
Niska R, Bhulya F, and Xu J. National Hospital Ambulatory Medical Care Survey: 2007 Emergency Department Summary. National health statistics reports; n°26. Hyattsville, MD: National Center for Health Statistics. 2010.

## **VIDAS®**, mini **VIDAS®** and **VIDAS® 3** A Solution adapted to Emergency Situations





# Emergency Diagnostics Solution

to Effective Triage of Emergency Department  
with Common Symptoms

of Breath - Fever - Cough\*



## CARDIAC

HF

NT-proBNP

VIDAS  
NT-proBNP

Galectin-3

VIDAS  
Galectin-3

ACS

Troponin I

VIDAS  
Troponin I Ultra

CK-MB

VIDAS CK-MB

Myoglobin

VIDAS Myoglobin

HF: Heart Failure - ACS: Acute Coronary Syndrome - PE: Pulmonary Embolism

### LARGE EMERGENCY TESTING PANEL AVAILABLE ON VIDAS® INSTRUMENTS\*

<b>CARDIAC</b>	VIDAS Troponin I Ultra	Ref. 30 448
	VIDAS CK-MB	Ref. 30 421
	VIDAS Myoglobin	Ref. 30 446
	VIDAS NT-proBNP	Ref. 30 449
	VIDAS NT-proBNP2	Ref. 30 458
	VIDAS Galectin-3	Ref. 41 1191
	VIDAS Digoxin	Ref. 30 603
<b>THROMBOSIS</b>	VIDAS D-Dimer Exclusion™ II	Ref. 30 455
<b>INFECTIONS</b>	VIDAS B.R.A.H.M.S PCT™	Ref. 30 450

\* Some of these reagents have not yet obtained regulatory clearance in some countries. Please contact your local bioMérieux representative for further information and product availability.



→ **Also available:**  
**Practical slide-rules for interpretation  
of NT-proBNP and D-Dimer cut-offs**  
**Clinical booklet : D-Dimer for exclusion of VTE**



→ **For further information, visit:**  
**[www.biomerieux.com/emergency](http://www.biomerieux.com/emergency)**

The information in this booklet is given as a guideline only  
and is not intended to be exhaustive.  
It in no way binds bioMérieux S.A. to the diagnosis established  
or the treatment prescribed by the physician.

bioMérieux S.A.  
69280 Marcy l'Etoile  
France  
Tel. : 33 (0)4 78 87 20 00  
Fax : 33 (0)4 78 87 20 90  
[www.biomerieux.com](http://www.biomerieux.com)

